

Biolumi[®] CX8

Integrated System



**Flexible & Compact System
Simplify Your Lab Operations**

Functional Modules



Multidimensional Indicator Lights

- Consumable indicators: confirming if the consumables remain sufficient
- Status indicator: Normal, Alarming or Error status will be shown in different colors




Built-in Operating Touchpad

- Real-time monitoring of reagents, consumables and waste
- Replacing reagents and consumables through the touch screen




Intelligent Software

- Smart maintenance with step-by-step guidance
- Bidirectional LIS connection by TCP/IP or COM
- Result system including functions of viewing, searching, deleting, detail viewing, retesting, uploading to LIS, recalculating and printing

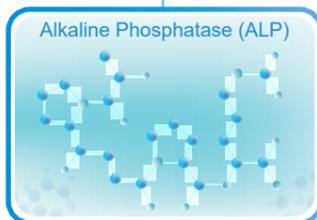

Biossays® C8 Module

- Throughput: 1600 T/H (BIO), 300 T/H (ISE)
- Bevel probe touches bottom to spit sample, ensuring the sampling accuracy of low-volume sample
- 10-step washing technology ensures thoroughly washing
- Full spiral mixing with adjustable speed
- Concave holographic grating post-splitting technology: high resolution and low stray light

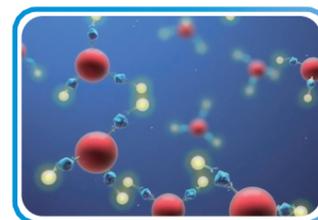


MAGLUMI® X8 Module

- Throughput: 600 T/H
- Up to 2800 walk-away tests
- No cross-contamination sampling and one pipetting for multiple tests technology
- 360° constant incubation technology
- 4-step high-efficiency washing unit
- Stable and accurate measurement



Decapper Module

- General throughput up to 400 T/H
- Safe operation with fully enclosed design
- UV irradiation in one-click, easy to operate
- Automatically decap and re-cap with the original tubes
- High throughput and space efficiency, suitable for different hospitals and laboratories

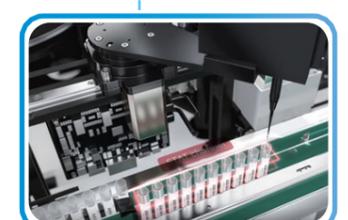

Sample Area

- 300 sample positions
- Supporting multiple sample types
- Storing and transporting samples
- Continuous loading, STAT position can be self-defined
- Indicator light shows sample status

Advantages of Biolumi® CX8



CONVENIENCE

- Intelligent consumable loading technology to reduce manual work
- Multiple built-in touchpads on each module offer easy operation
- Multidimensional indicator lights allow operators to check the status far away from the analyzer



RELIABILITY

- Proven magnetic microbeads separation and ABEI labeling techniques
- High-quality grating technology ensures stability and consistency
- Sample pipetting with disposal tips to avoid cross-contamination (Immunoassay)



EFFICIENCY

- Super high throughput to meet large-volume testing demands



Up to 600 tests/h
Immunoassay module



Up to 1600 tests/h
Biochemistry module

- Any sample position can be defined as STAT to process emergency samples efficiently
- Large number of reagent positions to deliver multiple tests in a single run



FLEXIBILITY

- Biolumi® CX8 Integrated System offers multiple configurations to meet diverse demands
- Decapper module configurable, no manual operation required
- Capable to link with Laboratory Automation system (TLA/LAS)



C8+X8			
	Biochemistry	Immunoassay	ISE
Throughput	1600 T/H	600 T/H	300 T/H
Dimension	326*118*150 cm		



C8*2 + X8*2			
	Biochemistry	Immunoassay	ISE
Throughput	3200 T/H	1200 T/H	600 T/H
Dimension	608*118*150 cm		



C8 + X8 + Decapper				
	Biochemistry	Immunoassay	ISE	Decapper
Throughput	1600 T/H	600 T/H	300 T/H	500 tubes/h (Decapping) 400 tubes/h (Re-capping)
Dimension	385*118*150 cm			



C8*2 + X8*2 + Decapper				
	Biochemistry	Immunoassay	ISE	Decapper
Throughput	3200 T/H	1200 T/H	600 T/H	500 tubes/h (Decapping) 400 tubes/h (Re-capping)
Dimension	666*118*150 cm			

One of the Broadest Test Menu

Immunoassay menu with 236 parameters:

Thyroid	Fertility	Autoimmune		
TSH (3rd Generation) ADV FT3 ADV TRAb T4 ADV Tg (Thyroglobulin) TMA T3 ADV TGA (Anti-Tg) Rev T3 FT4 ADV Anti-TPO T-Uptake	FSH free Testosterone SHBG LH DHEA-S Androstenedione HCG/ β -HCG ADV Progesterone PIGF PRL (Prolactin) free Estradiol sFlt-1 Estradiol 17-OH Progesterone *Inhibin A Testosterone AMH *Inhibin B	Connective Tissue Disease ANA Screen ENA Screen Anti-dsDNA IgG Anti-Sm IgG Anti-Rib-P IgG Anti-SS-B IgG Anti-SS-A IgG Anti-Jo-1 IgG Anti-Scl-70 IgG Anti-Centromeres IgG Anti-Histones IgG Anti-nRNP/Sm IgG *Anti-Ro-52 IgG *Anti-PM-Scl IgG *Anti-Nucleosome IgG	Rheumatoid Arthritis Anti-CCP *RF IgM *RF IgG *RF IgA *RF Screen	Endocrinology Anti-TPO TGA (Anti-Tg) TRAb TMA GAD 65 Anti-IA2 ICA IAA (Anti Insulin) *Anti-ZnT8
Bone Metabolism	Infectious Disease	Antiphospholipid Syndrome		
Calcitonin 25-OH Vitamin D ADV β -CTx Osteocalcin Intact PTH total P1NP *PTH (1-84)	Respiratory 2019-nCoV IgG 2019-nCoV IgM SARS-CoV-2 S-RBD IgG SARS-CoV-2 Neutralizing Antibody SARS-CoV-2 Ag <i>Mycoplasma Pneumoniae</i> IgG <i>Mycoplasma Pneumoniae</i> IgM <i>Chlamydia Pneumoniae</i> IgG <i>Chlamydia Pneumoniae</i> IgM *Respiratory Syncytial Virus IgM *Coxsackievirus B IgM *Adenovirus IgM *Influenza A Virus IgM *Influenza B Virus IgM *Legionella Pneumophila IgM *Human Parainfluenza Virus IgM	Autoimmune Liver Disease Anti-M2-3E IgG		
Hepatic Fibrosis	Cardiac	Inflammation Monitoring		
HA C IV Cholyglycine PIIIP N-P Laminin GP73	CK-MB NT-proBNP *sST2 Troponin I BNP *hs-cTnl STAT-EX™ Myoglobin D-Dimer *NT-proBNP STAT-EX™ hs-cTnl Lp-PLA2 *Myoglobin STAT-EX™ hs-CRP MPO *D-Dimer STAT-EX™ H-FABP HCY	CRP IL-6 (Interleukin 6) *PCT STAT-EX™ PCT (Procalcitonin) SAA (Serum Amyloid A) *CRP STAT-EX™ TNF- α *VEGF		
Immunoglobulins	Renal	Metabolism		
IgM IgA IgE IgG	Cr (CREA) Uric Acid (UA) Urea Cys-C α 1-MG β 2-MG mALB Urine / CSF Protein RBP NAG	Pepsinogen I Gastrin-17 IGF-I Pepsinogen II GH (hGH) IGFBP-3		
Glyco Metabolism	Cardiac	Coagulation Markers		
C-Peptide Anti-IA2 Proinsulin Insulin ICA *Glucagon GAD 65 IAA (Anti Insulin) *Anti-ZnT8	CK-MB NT-proBNP *sST2 Troponin I BNP *hs-cTnl STAT-EX™ Myoglobin D-Dimer *NT-proBNP STAT-EX™ hs-cTnl Lp-PLA2 *Myoglobin STAT-EX™ hs-CRP MPO *D-Dimer STAT-EX™ H-FABP HCY	D-Dimer TM tPAIC TAT PIC *FDP		
Tumor Markers	Cardiac	Prenatal Screening		
AFP PAP SCCA CEA CA 50 TPA-snibe Total PSA CYFRA 21-1 ProGRP f-PSA CA 242 HE4 CA 125 CA 72-4 HER-2 CA 15-3 NSE PIVKA-II CA 19-9 S-100 *AFP-L3%	CK-MB NT-proBNP *sST2 Troponin I BNP *hs-cTnl STAT-EX™ Myoglobin D-Dimer *NT-proBNP STAT-EX™ hs-cTnl Lp-PLA2 *Myoglobin STAT-EX™ hs-CRP MPO *D-Dimer STAT-EX™ H-FABP HCY	AFP (Prenatal Screening) free β -HCG *Inhibin A free Estriol PAPP-A		
Anemia	Cardiac	Hypertension		
Vitamin B12 EPO *Anti-Intrinsic Factor Ferritin RBC Folate *Active B12 Folate (FA)	CK-MB NT-proBNP *sST2 Troponin I BNP *hs-cTnl STAT-EX™ Myoglobin D-Dimer *NT-proBNP STAT-EX™ hs-cTnl Lp-PLA2 *Myoglobin STAT-EX™ hs-CRP MPO *D-Dimer STAT-EX™ H-FABP HCY	Direct Renin Angiotensin I Cortisol Aldosterone Angiotensin II ACTH		
STAT-EX™	Cardiac	Kidney Function		
*hs-cTnl STAT-EX™ *Myoglobin STAT-EX™ *PCT STAT-EX™ *NT-proBNP STAT-EX™ *D-Dimer STAT-EX™ *CRP STAT-EX™	CK-MB NT-proBNP *sST2 Troponin I BNP *hs-cTnl STAT-EX™ Myoglobin D-Dimer *NT-proBNP STAT-EX™ hs-cTnl Lp-PLA2 *Myoglobin STAT-EX™ hs-CRP MPO *D-Dimer STAT-EX™ H-FABP HCY	β_2 -MG Albumin *NGAL		
Drug Monitoring	Cardiac	Kidney Function		
Digoxin CSA (Cyclosporine A) FK 506 (Tacrolimus)	CK-MB NT-proBNP *sST2 Troponin I BNP *hs-cTnl STAT-EX™ Myoglobin D-Dimer *NT-proBNP STAT-EX™ hs-cTnl Lp-PLA2 *Myoglobin STAT-EX™ hs-CRP MPO *D-Dimer STAT-EX™ H-FABP HCY	β_2 -MG Albumin *NGAL		
Veterinary Testing	Cardiac	Kidney Function		
*cTSH *cTT4 *vFT4	CK-MB NT-proBNP *sST2 Troponin I BNP *hs-cTnl STAT-EX™ Myoglobin D-Dimer *NT-proBNP STAT-EX™ hs-cTnl Lp-PLA2 *Myoglobin STAT-EX™ hs-CRP MPO *D-Dimer STAT-EX™ H-FABP HCY	β_2 -MG Albumin *NGAL		

*Available soon

Biochemistry menu with 84 parameters:

Cardiac	Hepatic	Lipids	Inorganic Ion	Immune
CK CK-MB α -HBDH LDH Hcy LDH1 *D-Dimer *ACE *IMA *MYO	TBA ALT (SGPT) AST (SGOT) ALP GGT TBIL (Vanadate) DBIL (Vanadate) TP ALB Ammonia (AMM) PA CHE AFU 5'-NT Haptoglobin (HPT) *TBIL (DPD) *DBIL (DCA) *CG (Latex) *CG (CEDIA) *GLDH *ADA *CER *LAP *GR	HDL-C LDL-C TC TG ApoE ApoA1 ApoB Lp (a) NEFA	Ca (OCPC) P (PHOS) Mg *Ca (Arsenazo III) *CO2	IgM IgA IgG *C3 *C4
Diabetes	Renal	Renal	Anemia	ISE
GLU LAC HbA1c GSP D3H GA	Cr (CREA) Uric Acid (UA) Urea Cys-C α 1-MG β 2-MG mALB Urine / CSF Protein RBP NAG	Cr (CREA) Uric Acid (UA) Urea Cys-C α 1-MG β 2-MG mALB Urine / CSF Protein RBP NAG	Transferrin (TRF) Iron UIBC *G6PD *Ferritin	Na ⁺ K ⁺ Cl ⁻ Ca ²⁺ pH
Pancreatic	Renal	Renal	Rheumatism	Inflammation
α -AMY LIP	Cr (CREA) Uric Acid (UA) Urea Cys-C α 1-MG β 2-MG mALB Urine / CSF Protein RBP NAG	Cr (CREA) Uric Acid (UA) Urea Cys-C α 1-MG β 2-MG mALB Urine / CSF Protein RBP NAG	ASO RF	CRP (Full Range)

*Available soon

Technical Specifications

Sample Area	Sample type	Serum, Plasma, Whole blood, Urine, CSF
	Sample position	300
	Sample recognition	Barcode
Biochemistry Module	Throughput	1600 T/H
	Reaction cuvettes	362 cuvettes in total: 181 cuvettes each for two segments
	Reagent position	72 positions each for R1 and R2
	Sample volume	1.5 μ L ~ 25.0 μ L (0.1 μ L increments)
	Sample probe	Clot detection, liquid detection, collision detection and auto-clearing
	Reagent volume	20 μ L-250 μ L (1 μ L increments)
	Reagent probe	Liquid detection, collision detection and auto-clearing
	Reaction volume	80 μ L ~ 250 μ L
	Reaction time	10 min, 22 min
	Reaction temperature	37.0 $^{\circ}$ C \pm 0.2 $^{\circ}$ C (with the fluctuation not greater than 0.1 $^{\circ}$ C)
	Light Source	12 V, 100 W halogen lamp
	Wavelength accuracy	13 fixed wavelengths (340,380, 405, 450, 480, 505, 546, 570,600, 660, 700, 750, 800 nm)
Immunoassay Module	Throughput	Up to 600 T/H
	Walk-away tests/time	2800 Tests
	Reagent position	42
	Sample volume	10 μ L ~ 280 μ L
	Sample probe	Liquid detection, clot detection, bubble detection, collision detection and volume tracking
	Reagent volume	10 μ L ~ 450 μ L
	Reagent probe	Auto-clearing, Liquid detection, collision detection and volume tracking
	Reagent area	Working temperature: 8 $^{\circ}$ C ~ 12 $^{\circ}$ C Storage temperature: 2 $^{\circ}$ C ~ 8 $^{\circ}$ C
	Reaction cuvette number	2912/5824
	Reaction temperature	37.0 $^{\circ}$ C \pm 0.3 $^{\circ}$ C (with the fluctuation not greater than 0.1 $^{\circ}$ C)
	Disposable tip number	576/1152
ISE Module	Sample volume	90 μ L
	Information recognition	RFID
	Reagent type	Calibrator A, Calibrator B
	Electrode type	K ⁺ , Na ⁺ , Cl ⁻ , iCa ²⁺ , pH, reference electrode

MAGLUMI®, Biossays® and Biolumi® are trademarks of Snibe.

All other product names and trademarks are the property of their respective owners.



Shenzhen New Industries Biomedical Engineering Co., Ltd.
(SNIBE Co.,Ltd.)

No.23, Jinxiu East Road, Pingshan District, 518122 Shenzhen,
P.R. China

Tel: +86 755 26501514

Email: sales@snibe.com

Fax: +86 755 26654850

Web: www.snibe.com

Follow us on

